


D1.5 Oil-Free Compressor



100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
SKU: 2880116

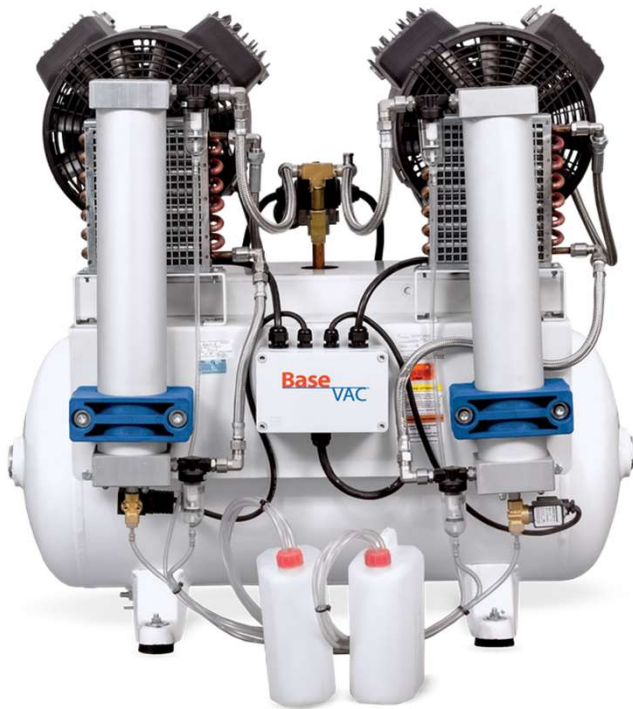

Part #	2880116
Max (HV) Users	6
Operatories	5-7
Hp	2 x 1.5
kW	2.3
Duty Cycle	60%
Amps	17
Volts	220V
Sound Level (dB)	74 dB
L x W x H	35" x 24" x 36"
Ibs.	231
Tank	20 Gal, ASME & CRN Rated
CFM Max	10 @ 80PSI
Operating PSI	80 - 100PSI
Dryer Type	Desiccant
Silent Options (-20db)	No

Product Highlights:



100% Clean, Dry Medical-Grade Air
100% oil-free with Perma Dry desiccant drying system.



Low Maintenance Cost
Simple Installation & Easy to Clean Inlet Filter.



Robust & Powerful Motors
Long lasting Italian made motors



2-Year Warranty
Unlimited hour hassle-free warranty.

Package Includes:

- Dual Motor Compressor
- ASME & CRN Rated 20-Gal Tank
- Dual Desiccant Dryers with Perma Dry Tech
- Pressure Filter Regulator
- Auto Condensate Drain
- 6' length of 1/4" Flex Line hose


Optional Accessories:

Part #	Description
2900601	Buck/Boost Transformer
2900107	Remote Start Transformer
2880131	External Air Inlet


1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com

D1.5 Oil-Free Compressor

100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
SKU: 2880116

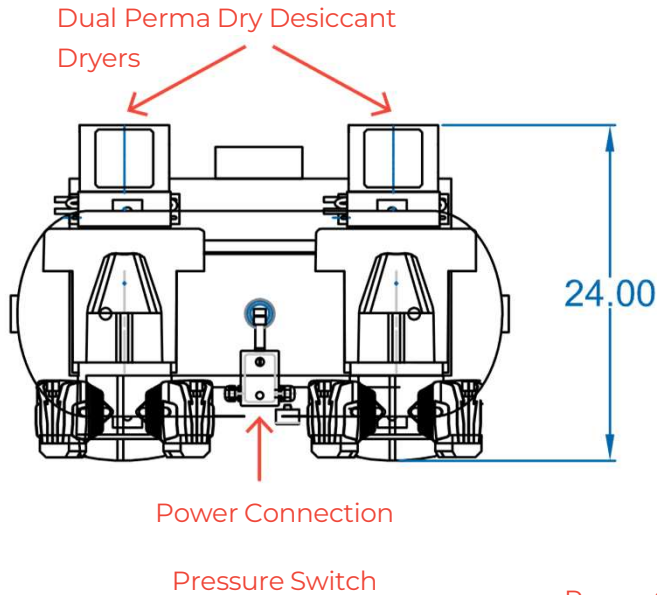

Site Requirements	
Voltage (Min/Max)	208 / 240V 1 60Hz
Full Load Amps	17
Breaker Rating	30 Amps
Wire Size	10
Watts/Hour	1286
BTU/Hour	4500
Install Footprint	47" x 36" x 41"
Shipping Weight	286 lbs.

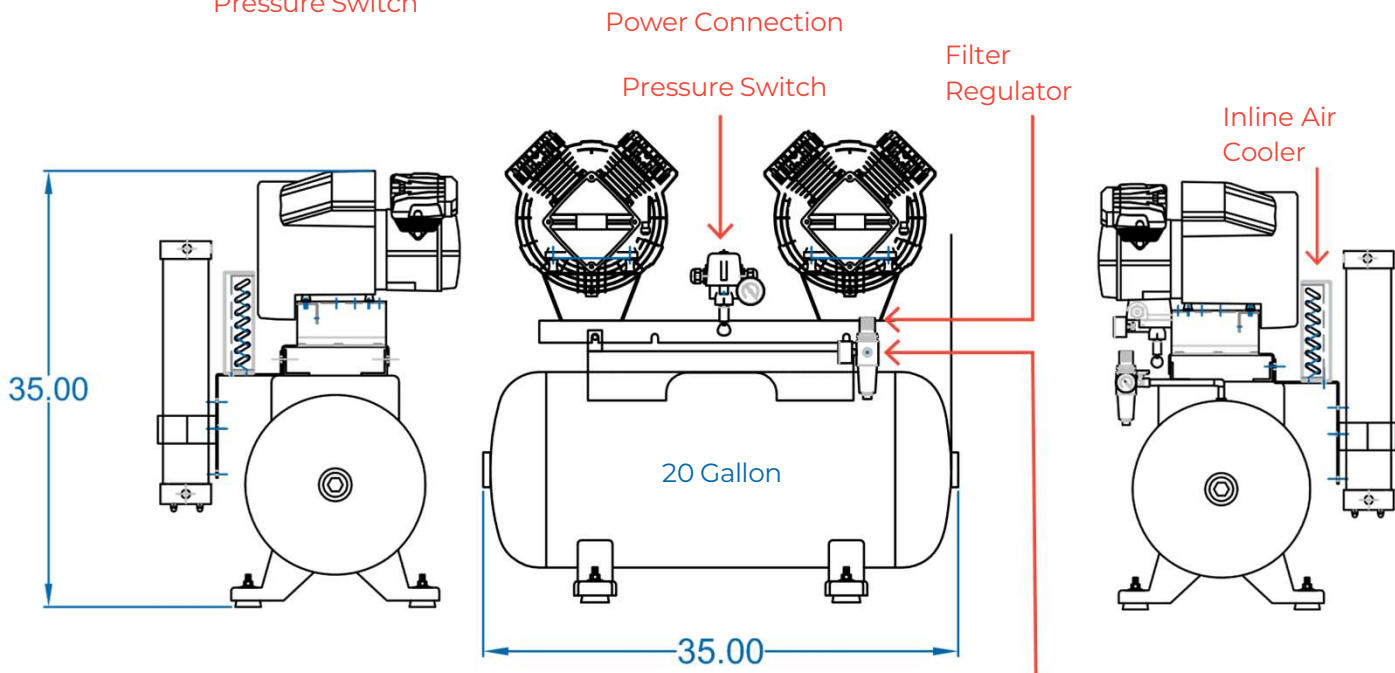

Note:

- Clean air intake manifold combines compressor heads to a common 3/4" female copper fitting. For clean air intake, use 1.5" - 2" PVC or ABS schedule 40. Hospital and government clinics may require copper / stainless steel pipe (see NFPA or CSA requirements).
- Recommended to allow 18" around and above the compressor for easy access service
- Room temperature should range 50 °F - 104 °F (10 °C - 40 °C).

Hose Connection

Hose Connection to main at regulator - 3/8 FNPT - Compressor supplied with 6-foot flexible pressure hose for connection to 1/2-inch type "L" of Type "K" copper pipe in accordance with medical air requirements - refer to NFPA or CSA for more detail.

1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com