


D2.0 Oil-Free Compressor



100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
SKU: 2880111

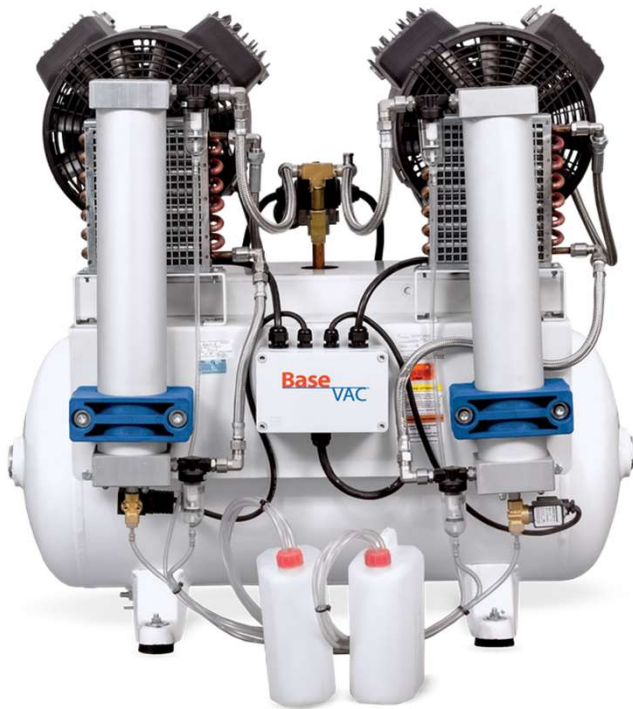

Part #	2880111
Max (HV) Users	8
Operatories	8-10
Hp	2 x 2
kW	3.3
Duty Cycle	60%
Amps	18
Volts	220V
Sound Level (dB)	74 dB
L x W x H	35" x 24" x 36"
Ibs.	242
Tank	20 Gal, ASME & CRN Rated
CFM Max	14 @ 80PSI
Operating PSI	80 - 100PSI
Dryer Type	Desiccant
Silent Options (-20db)	No

Product Highlights:



100% Clean, Dry Medical-Grade Air
100% oil-free with Perma Dry desiccant drying system.



Low Maintenance Cost
Simple Installation & Easy to Clean Inlet Filter.



Robust & Powerful Motors
Long lasting Italian made motors



2-Year Warranty
Unlimited hour hassle-free warranty.

Package Includes:

- Dual Motor Compressor
- ASME & CRN Rated 20-Gal Tank
- Dual Desiccant Dryers with Perma Dry Tech
- Pressure Filter Regulator
- Auto Condensate Drain
- 6' length of 1/4" Flex Line hose


Optional Accessories:

Part #	Description
2900601	Buck/Boost Transformer
2900107	Remote Start Transformer
2880131	External Air Inlet


1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com

D2.0 Oil-Free Compressor

100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
SKU: 2880111

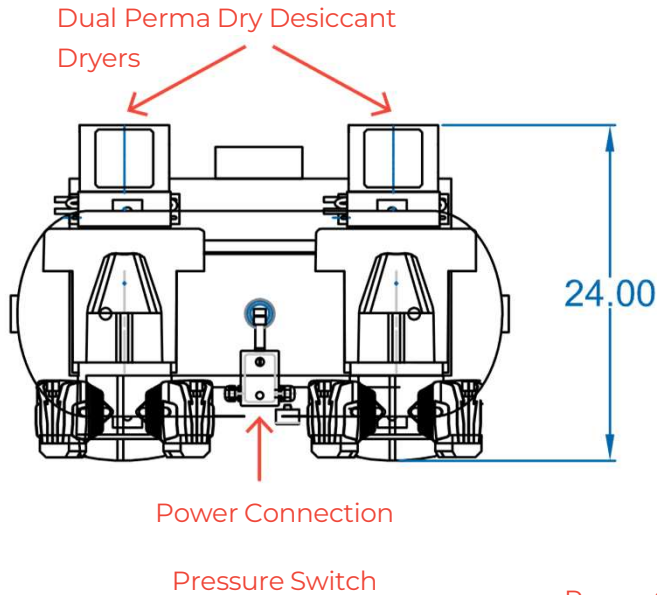

Site Requirements	
Voltage (Min/Max)	208 / 240V 1 60Hz
Full Load Amps	18
Breaker Rating	30 Amps
Wire Size	10
Watts/Hour	1485
BTU/Hour	6000
Install Footprint	47" x 36" x 41"
Shipping Weight	297 lbs.

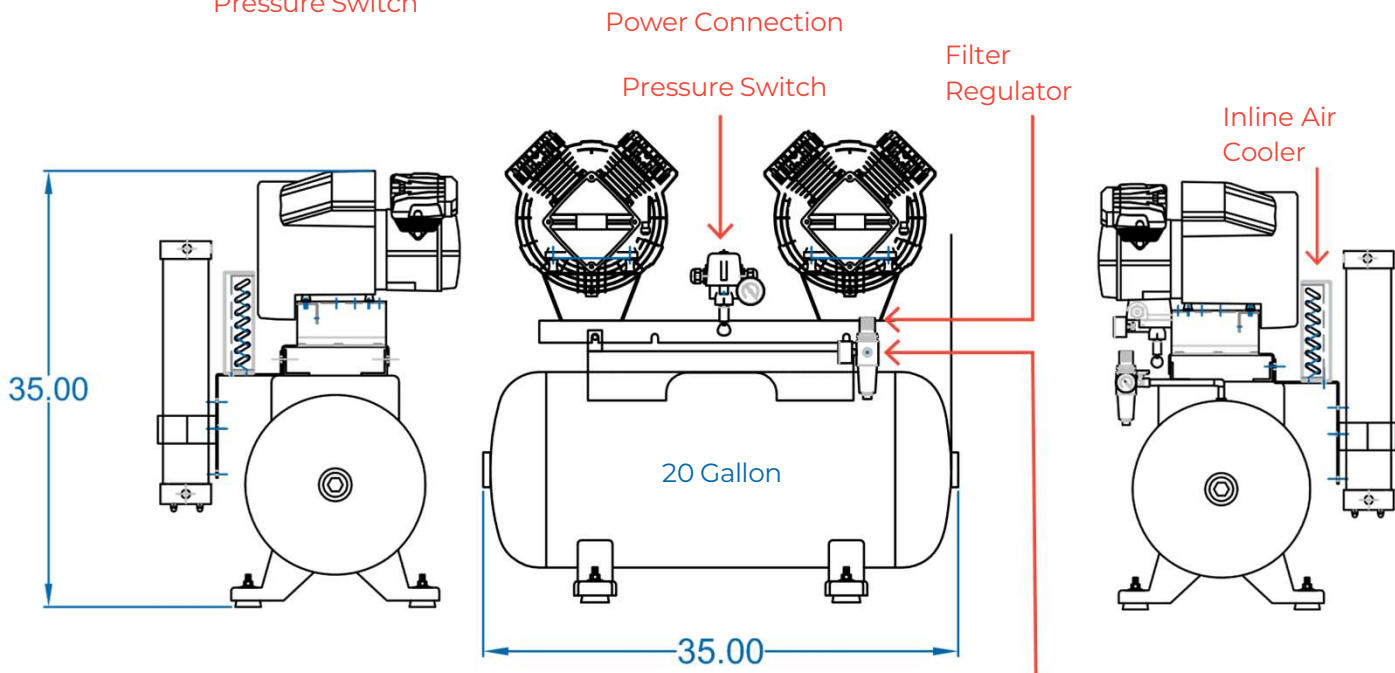

Note:

- Clean air intake manifold combines compressor heads to a common 3/4" female copper fitting. For clean air intake, use 1.5" - 2" PVC or ABS schedule 40. Hospital and government clinics may require copper / stainless steel pipe (see NFPA or CSA requirements).
- Recommended to allow 18" around and above the compressor for easy access service
- Room temperature should range 50 °F - 104 °F (10 °C - 40 °C).

Hose Connection

Hose Connection to main at regulator - 3/8 FNPT - Compressor supplied with 6-foot flexible pressure hose for connection to 1/2-inch type "L" of Type "K" copper pipe in accordance with medical air requirements - refer to NFPA or CSA for more detail.

1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com