

S2.0 Oil-Free Compressor



100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
SKU: 2880109



Part #	2880109
Max (HV) Users	4
Operatories	3-5
Hp	2
kW	1.5
Amps	9
Volts (Min-Max)	220V
Sound Level (dB)	74 dB
Tank LxWxH	26" x 25" x 30"
Tank Weight	154
Tank	13 Gal, ASME & CRN Rated
CFM Max	6.5 @ 80PSI
Operating PSI	80-100 PSI
Dryer Type	Desiccant
Silent Options (-20db)	Yes



Package Includes:

- Compressor
- Tank
- Perma Dry Desiccant Dryer
- Pressure Filter Regulator
- Auto Condensate Drain
- 6' length of 1/4" Flex Line hose

Optional Parts:

- buck Boost Part No: 2900601
- Remote Start Transformer Part No: 2900107
- External Air Inlet Part No: 2880131

Note:



- Clean air intake manifold combines both compressor heads to a common 3/4" female copper fitting. For clean air intake from fresh air source use 1.5" - 2" PVC or ABS schedule 40 or better. Hospital and government clinics may require copper / stainless steel pipe (see NFPA or CSA requirements).
- Room must be provided with active ventilation. It's recommended to allow 18" around and above the compressor for easy access service.
- Room temperature should range 50 °F - 104 °F (10 °C - 40 °C).

1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com

S2.0 Oil-Free Compressor


100% Oil Free | Dry Medical Air | Unlimited Hour Warranty
 SKU: 2880109

Site Requirements	
Voltage (Min/Max)	208 / 240V 1 60Hz
Full Load Amps	9
Breaker Rating	20 Amps
Wire Size	12
Watts/Hour	825
BTU/Hour	3000
Install Footprint	39" x 35" x 38"
Shipping Weight	219 lbs.

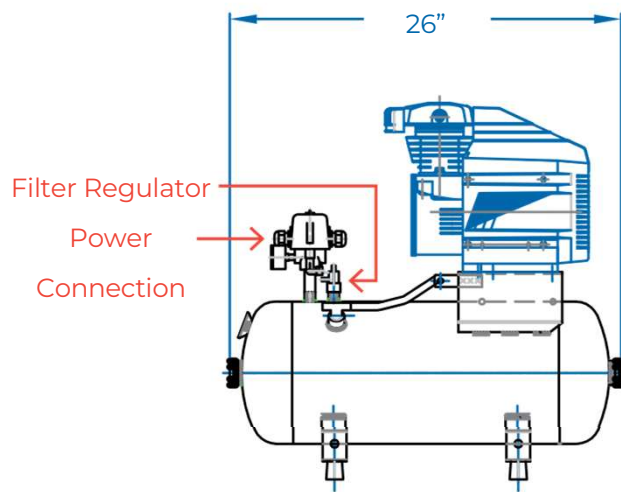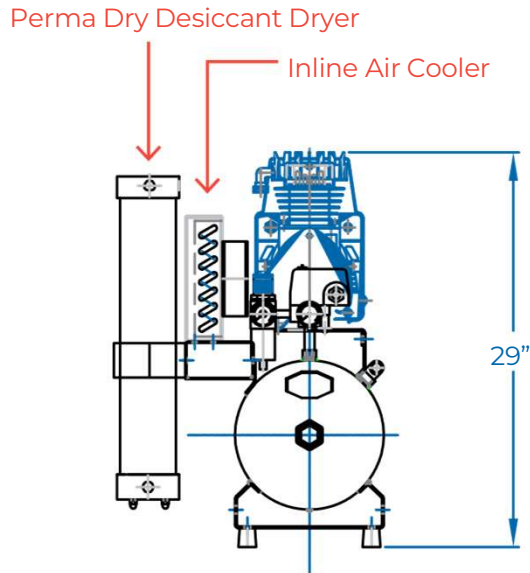

Hose Connection

Hose Connection to main at regulator - 1/4" FNPT - Compressor supplied with 6 foot flexible pressure hose for connection to 1/2 inch type "L" of Type "K" copper pipe in accordance with medical air requirements - refer to NFPA or CSA for more detail.

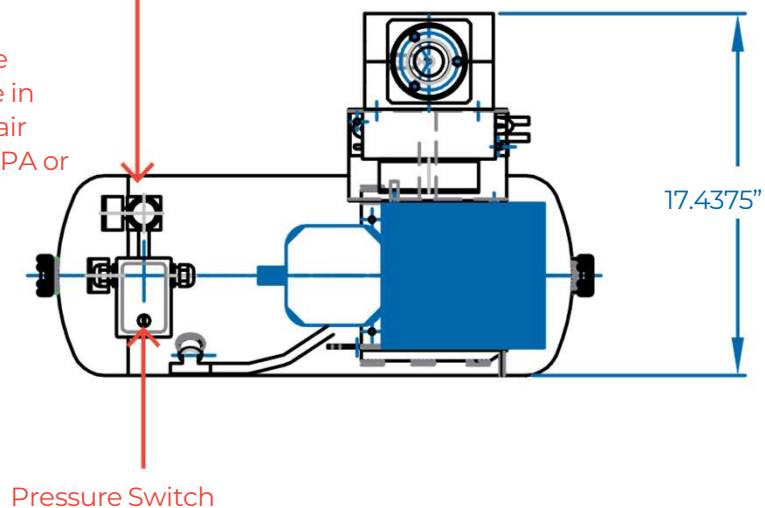

1.800.668.8736 - sales@basevacdental.com - www.basevacdental.com